UMAX U-SONIC WHITE

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» CLEAN

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Speed 34 800 strokes/min.

» WHITE

Mode for teeth whitening. First two minutes for plaque removal followed by 30 seconds for teeth polishing.

Speed 33 000 - 40 000 strokes/min.

» POLISH

Mode for a shining smile and teeth polishing.

Speed 35 000 strokes/min.

» SOFT

Soft cleaning of teeth and gums.

Speed 33 000 strokes/min.

» SENSITIVE

Mode for cleaning of sensitive teeth and gums.

Also good for beginners.

Speed 28 200 strokes/min.

••• HOW TO USE

» 1

Connect the head of the toothbrush to the body.

» 2

Wet the toothbrush bristles and put on a toothpaste.

» 3

Put the toothbrush next to your teeth and then turn it on.

» 4

Slowly move the toothbrush head over all teeth. It is not necessary to move the toothbrush up and down and left and right like you would do using a regular toothbrush. Every 30 seconds there is a short pause to indicate that you should change position.

••• CLEANING OF TOOTHBRUSH

» 1

Clean your toothbrush regularly. Wash the toothbrush head under running water.

» 2

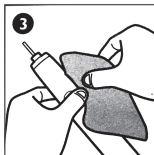
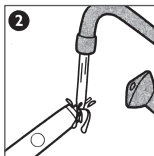
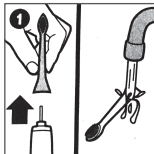
Wash the toothbrush body in the same way.

» 3

Dry the toothbrush head, body, and charging base using a soft cloth. You can also use soap.

» 4

Do not put the toothbrush body and charging base into a dishwasher.

••• CHARGING

» 1

Charge the toothbrush for 12 hours before first use.

» 2

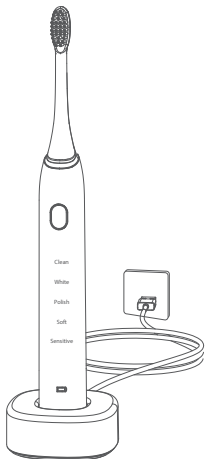
The toothbrush will vibrate to indicate low battery power. The LED battery indicator will lit up.

» 3

Put the toothbrush into the charging base. The LED indicator will be blinking during charging.

» 4

When the battery is fully charged, the battery indicator will lit up.

••• WARNINGS

» 1

Keep charging base away from water. Do not place it near wet areas, such a washbasin. Do not put the charging base under water. After cleaning, make sure that the charging base is completely dry before plugging it in.

» 2

Do not put the toothbrush head, body, and charging base into dishwasher.

» 3

Do not use the toothbrush in a shower or bathtub.

» 4

If you have had a teeth or gum operation in the previous two months, consult your dentist first about using the toothbrush.

» 5

Replace the toothbrush head every three months.

••• SUPPORT**» TECHNICAL SUPPORT**

Email: support@umax.cz

» SERVICE CENTER

ConQuest entertainment a.s.

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198 00 Praha 9

+420 284 000 111

Opening Hours:

mon-thu 9:00 - 17:00

fri 9:00 - 16:00

UMAX

© 2018 Umax Czech a.s.